# FLAEM<sup>®</sup>

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## Aerosol therapy unit

*Congratulations on your purchase and thank you for choosing our unit. Our goal is to fully satisfy consumers by offering them cutting-edge systems for treating ailments of the upper airways.*

***Carefully read these instructions and keep them in a safe place for future reference.***

***Use the unit only as described in this instruction manual.***

***This is a medical device for home use and may only be used with medication prescribed or recommended by a doctor.***

*Visit our internet site to view the whole range of Flaem products: [www.fluemnuova.it](http://www.fluemnuova.it)*

## INSTRUCTIONS FOR USE

**Before each use, the nebuliser and accessories should be cleaned according to the instructions in the section "Cleaning, disinfection and sterilisation". A personal use of the nebuliser and its accessories is recommended in order to avoid any risks of infection.**

1. Plug the power cord into a power socket corresponding to the voltage of the unit.

- ⚠ • Before using the unit for the first time, and periodically throughout its life, check the power supply cord for any damage; if you detect any damage, do not plug it in. Take the appliance immediately to an authorised service centre or your dealer.
- Never use adapters for power supplies different from the one indicated on the bottom of the unit (plate details). Keep the cord away from heated surfaces.

2. Carefully wash your hands before preparing the inhalation therapy.

- ⚠ Do not handle the plug with wet hands and do not use the appliance while bathing or showering.

3. Open the nebuliser by turning the upper part (9c.) counterclockwise.

4. Pour the medication prescribed by the physician into the lower part (9a.). Close the nebuliser by turning the upper part (9c.) clockwise.

5. Connect the accessories as shown in the "Assembly diagram". Sit comfortably, holding the nebuliser in your hands, place the mouthpiece over your mouth or use the nasal prong or mask.

- ⚠ **Interactions:** the materials used that come into contact with drugs have been tested with a wide range of drugs. However, due to the variety and continuous evolution of medication, the possibility of a chemical interaction cannot be excluded. Once opened, the medication should be used as soon as possible, and prolonged contact between drug and nebuliser should be avoided.

6. Start the unit by turning the switch on (.6.). Cover the hole of the nebulisation manual control (.8.) with your finger and breathe smoothly and deeply; after inspiration we recommend to hold your breath for a moment in order to allow inhaled aerosol droplets to deposit; whilst holding your breath remove the finger from the hole of the nebulisation manual control (.8.) to optimise drug intake. Breathe out slowly.

**To obtain continuous nebulisation, do not fit the nebulisation manual control (.8.), especially in the case of children or invalids.**

- ⚠ • Children and disabled persons should always use the unit under strict supervision of an adult who has read this manual.
- The unit should not be used in the presence of anaesthetic mixture inflammable with air, oxygen or nitrous oxide.
- Do not obstruct the air slots placed on both sides of the unit. During use, always place the unit on a firm surface clear of obstacles. The air slots should be kept free of obstruction by gauze, hairs or similar objects. Do not place any objects inside the air slots.

7. When the treatment is finished, turn off and unplug the unit.

- ⚠ • Always unplug the power cord immediately after use.

8. It may happen that during the use of the nebuliser visible humidity deposits form in the connection tube (.7.); disconnect the tube from the nebuliser and let it dry using the compressor air flow in order to reduce any risks of infection between users.

## CLEANING, DISINFECTION and STERILISATION

*Before undertaking any cleaning operation, switch off and unplug the unit.*

## Accessories

- **Cleaning:** Wash under running water with mild (non abrasive) washing-up liquid or in the dishwasher on warm cycle, after disassembling all the parts. Clean the nozzle (9b.) after detaching it from the upper part (9c.) of the nebuliser and disassembling as shown in fig. A.
- **Disinfection:** All accessories (except for the connecting tube) can be disinfected choosing one of these methods: a) by soaking in a solution of water (60%) and white vinegar (40%). Finally, rinse thoroughly with drinking water; b) by deep boiling (20 minutes); this method cannot be used for the tube and the mask.
- **Sterilisation:** Sterilise all parts with cold chemical sterilizers (peracetic acid based, following the instructions provided by the manufacturer), or (except for the tube and the mask) in an autoclave at 121°C (1,2 bar) for 20 minutes or at 134°C (1.2 bar) for 8 minutes max. carefully following the instructions provided by the autoclave manufacturer.
- **Drying:** After each use dry all parts very carefully with a paper tissue or a warm air flow (hair dryer) in a non-humid place. Reassemble the nozzle and fit it back in place on the nebuliser as shown in the "Assembly diagram".

⚠ When used very frequently, we recommend that the nebuliser be replaced every 6 months (or earlier if it becomes clogged), in order to guarantee the maximum therapeutic effect.


### Unit and exterior of the tube

To clean use only a damp cloth and an antibacterial detergent (non-abrasive and solvent free).

- ⚠
- The casing of this unit is not protected against the penetration of liquids. Never wash the appliance under running water or by soaking.
  - Never submerge the unit in water; should this accidentally occur, unplug immediately. Do not remove or touch the appliance under water before pulling out the plug. Take it immediately to an authorised service centre or your dealer.


## RF6 PLUS NEBULISER DESCRIPTION

It is professional, quick and suitable for administering all types of drugs, including more costly ones, also in patients with chronic diseases. It features a **valve system** that is automatically controlled by the patient's breathing. The configuration of the internal lines of the RF6 PLUS nebuliser ensures the ideal granulometry providing effective treatment right down to the lower airways.

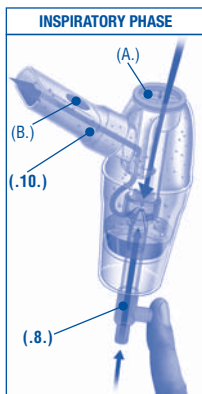During the inspiratory phase, placing a finger over the hole of the Nebulisation manual control causes compressed air to be introduced into the nebuliser; at the same time the inspiratory valve (A.) of the nebuliser opens to permit entry of an additional flow of air. The combination of these air flows increases the quantity of medication nebulised and thus assimilated by the patient (MAXIMUM ASSIMILATION OF MEDICATION). During the whole inspiratory phase the expiratory valve (B.) remains closed.

In the expiratory phase, removing the finger from the Nebulisation Manual Control will cut off the flow of compressed air into the nebuliser and thus interrupt the nebulisation of medication; at the same time the inspiratory valve (A.) will close, preventing the nebulised mist from being discharged from the nebuliser and dispersed in the surrounding environment (MINIMUM WASTE OF MEDICATION). The air exhaled by the patient will be released through the valve (B.), which opens automatically during the expiratory phase while the mouthpiece remains in the mouth.
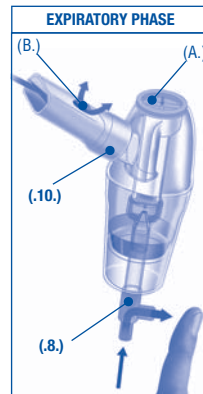
### NON-SPILL SYSTEM

The RF6 PLUS nebuliser is specially suitable also for pediatric use, since the non-spill system prevents medication from spilling out even if the nebuliser is accidentally turned upside-down.

(Efficacy tested up to 5 ml, the average dose recommended for normal use)


## AIR FILTERING

The unit is equipped with an air filter (.4.) which should be replaced when dirty or when its colour changes. Do not wash or re-use the used filter. The filter must be regularly replaced to help ensure efficient compressor performance. The filter should be regularly checked. For spare filters contact your dealer or an authorised service centre.

### To replace the filter:

Insert a flat screwdriver between the filter edge and the unit body. Lift the filter and remove it by rotating and pulling it upwards. The filter is designed to remain fixed in its housing.

**Use only Flaem's original accessories.**
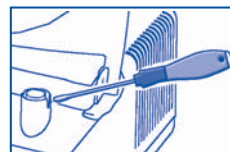
## SPARE PARTS

### Description

- RF6 PLUS nebuliser kit
- 1 m Connection tube
- PVC Adult mask
- PVC Child mask
- Spare Air filter kit per unit (n°2 pcs)

### Code

AC0234P  
AC035  
AC079  
AC089  
AC0164P

Elisir can be used with **RHINO CLEAR** Nasal wash

## TECHNICAL FEATURES

### Aerosol therapy unit Mod. ELISIR F1000

Voltage:	230V ~ 50Hz 140VA	115V ~ 60Hz 180VA	220V ~ 60Hz 160VA	100V ~ 50/60Hz
Fuse:	T 2A - 250V	T 3A - 250V	T 2A - 250V	T 5A - 250V
In compliance with:	Dir. 93/42/EEC (EN 60 601-1/A2) <b>CE</b> 0051			
Approvals:	<b>KEMA</b> EUR			

Max pressure: 2.5 bar approx.  
Compressor air output: 10 LPM approx.  
Dimensions: 20x30x10 cm  
Weight: 2.100 Kg  
Sound level (at 1 m): 57 dB (A)  
Continuous use

### RF6 PLUS nebuliser

Medication minimum capacity: 2 ml  
Medication maximum capacity: 8 ml  
Operating pressure (with neb.): 0.9 bar approx.  
Max delivery<sup>(1)</sup>: 0.42 ml/min approx.  
MMAD<sup>(1) (2)</sup>: 1.53 µm  
Particles % < 5 µm<sup>(1) (2)</sup>: 99 %

<sup>(1)</sup> Average values measured with 0.9% saline solution in compliance with Flaem internal procedure.  
<sup>(2)</sup> Certified with TSI Aerosizer 3220 computerised laser system in compliance with Flaem internal procedure.

**In vitro characterization carried out at the Inamed Research GmbH & Co. KG on behalf of TÜV Rheinland Product Safety GmbH – Germany in compliance with the new European Standard for aerosol therapy appliances, Standard EN 13544-1. Further details are available upon request.**

- MMAD: 3,2 µm  
- Breathable fraction < 5 µm (FPF): 73 %

**Operating requirements:** temperature min. 10°C; max 40°C  
**RH** air humidity min. 10%; max 95%

**Storage requirements:**  temperature min. -25°C; max 70°C  
**RH** air humidity min. 10%; max 95%

**Operating/storage atmospheric pressure:** **P** min. 690 hPa; max. 1060 hPa

-  • The Manufacturer, the Retailer and the Importer are responsible for safety, reliability and performance only if: a) the unit is used following the instructions for use; b) the electrical system of the premises where the appliance is used is in conformity with current laws.  
• Repairs should be performed by authorised personnel only. Unauthorised repairs will invalidate the warranty and may represent a hazard to the user.

 **Disposal:** In conformity with Directive 2002/96/EC on waste electrical and electronic equipment (WEEE), at the end of its life the unit should be taken to an electrical appliance collection point for its environmentally friendly disposal.

## SYMBOLS

 Class II unit

 B-type unit

 Important safeguards  
Warning: follow the operating instructions

 Alternating current

### Electromagnetic compatibility:

**ELISIR** unit is designed to meet the current requirements for electromagnetic compatibility. However if you suspect that the unit is interfering with the operation of your TV, radio or other electrical unit, try relocating the unit until the interference ceases. Alternatively try plugging the unit into a different power socket.

The Manufacturer reserves the right to make technical and functional changes to this product without prior notice.

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